

ZEISS MultiSEM

The World's Fastest Scanning Electron Microscope

zeiss.com/multisem



Seeing beyond


Revolutionize the Speed of Electron Microscopy with ZEISS

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › Technology and Details

Unleash the speed of this unique scanning electron microscope and start thinking about new dimensions. Now, at last, you can image huge samples at nanometer resolution, driven by the unrivalled acquisition speed of MultiSEM.

MultiSEM is designed for continuous, reliable 24/7 operation. Simply set up your high-throughput data acquisition workflow. Then get on with your day while MultiSEM takes high contrast images all by itself – no supervision needed.

MultiSEM uses ZEN imaging software, so you can control this powerful microscope in an intuitive yet flexible way. Automated tuning routines make sure you achieve the best high resolution data.

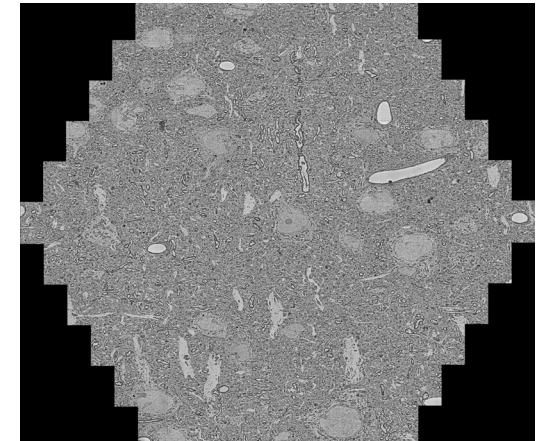

Mouse brain section, image acquired using ZEISS MultiSEM 506 with 91 beams and at 4 nm pixel size. The hexagonal field of view (FoV) is 165 μm wide. Sample: courtesy of J. Lichtman, Harvard University, Cambridge, MA, USA

Leading Edge Science

ZEISS offers you access to a groundbreaking technology for a rapidly developing application space. Take profit from the opportunities of a new technology and become member of a worldwide growing community at the forefront of science.

By providing individual application support and a Premium Service concept, we help you to explore new territory. Are you ready to take the next step?

Simpler. More Intelligent. More Integrated.

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details

Highest-ever Acquisition Speed at Nanometer Resolution


Multiple electron beams working in parallel give you unprecedented gross imaging speed. Acquiring an area of 1 mm^2 at 4 nm pixel size takes only a few minutes of imaging time. The unrivaled imaging speed of more than 1 TB per hour enables imaging of large volumes ($> 1 \text{ mm}^3$) at nanometer resolution. Optimized detectors collect the secondary electron signals very efficiently, providing you with high contrast images at low noise levels.

Electron Microscopy of Huge Samples


MultiSEM is built for continuous 24/7 operation and equipped with a sample holder covering an area of $10 \text{ cm} \times 10 \text{ cm}$. That means you no longer have to sacrifice sample size for nanometer resolution. You can finally image the entire sample and discover everything you need to answer your scientific questions. With automated acquisition protocols to enable large area imaging, you will get the detailed full picture, without losing the macroscopic context.

Electron Microscopy with ZEN Imaging Software


By introducing ZEN to MultiSEM, we bring the standard software for ZEISS light microscopes to the world of electron microscopy. ZEN lets you control MultiSEM in a straightforward, intuitive way. Smart auto-tuning routines support you as you capture optimal images with high resolution and quality. You quickly and easily set up even complex automated acquisition procedures, adapted and tuned to your sample imaging. ZEN for MultiSEM also masters the high speed required for continuous, parallel image recording. An open application programming interface (API) is provided for flexible and fast application development.

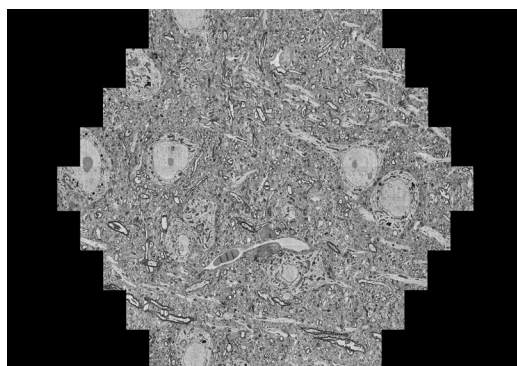

Mouse brain section (50 nm thick) image acquired in 1.3 s with a MultiSEM covering a hexagonal field of view of $108 \mu\text{m} \times 94 \mu\text{m}$. Sample: courtesy of J. Lichtman, Harvard University, Cambridge, MA, USA

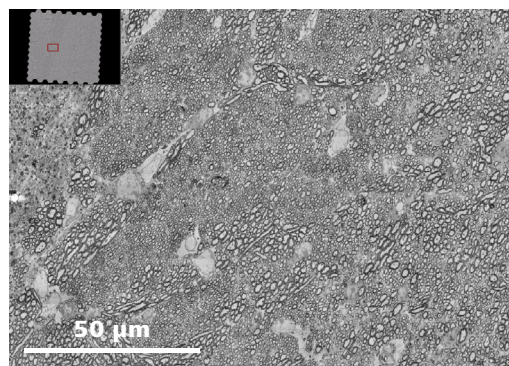


Mouse brain section, automated image acquisition of 1 mm^2 at 4 nm pixel size totaling 100 GB of data. Sample: courtesy of J. Lichtman, Harvard University, Cambridge, MA, USA

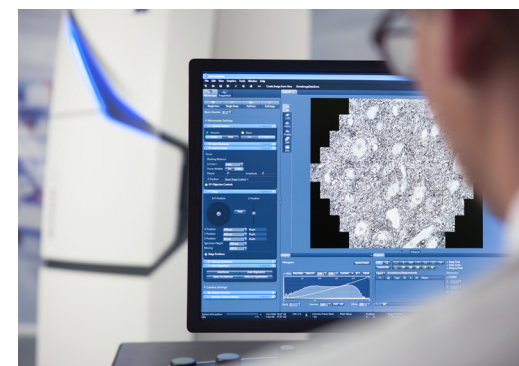


ZEN imaging software is used to control the MultiSEM.

Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details

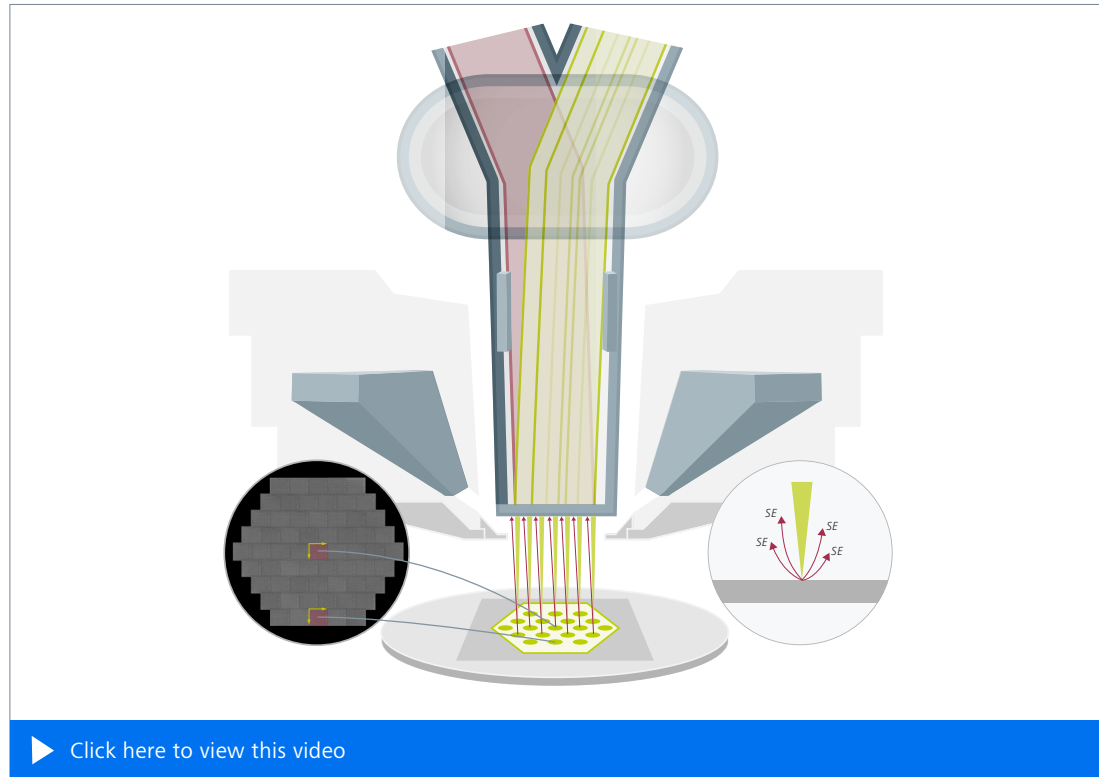

MultiSEM achieves high imaging speed by employing multiple electron beams and detectors in parallel. The key to this approach is a finely tuned detection path (red) collecting a large yield of secondary electrons used for imaging with a multiple detector array. Each beam carries out a synchronized scanning routine at one sample position, resulting in a single sub-image. The electron beams are arranged in a well-characterized hexagonal pattern. By merging all sub-images together, the final, full image is formed. A parallel computer setup is used for fast data recording ensuring high total imaging speed. Image acquisition and workflow control are fully separated in the MultiSEM system to guarantee full performance.

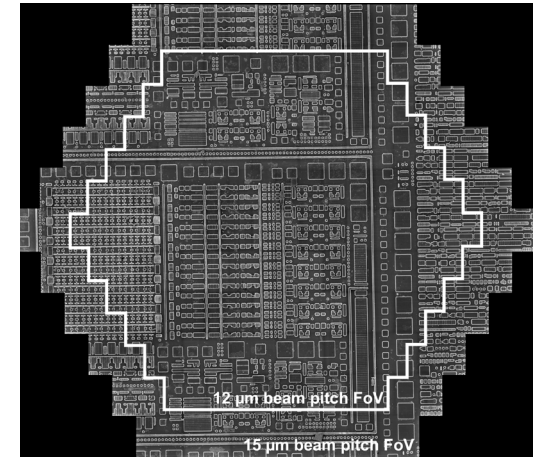

MultiSEM 506 operates 91 electron beams in parallel, covering a multi-beam field of view of at least 132 μm and a maximum of 165 μm width. For more throughput, the system can optionally be configured with a wider beam-to-beam distance. You can choose between the standard 12 μm beam pitch configuration and the extended 15 μm beam pitch configuration.

The image above (graphics processor chip) shows the fields of view of the two beam pitch options available in comparison. With every single scan pass, MultiSEM 506 configured with a beam-to-beam distance of 15 μm covers an area more than 50% larger as compared to the 12 μm beam-to-beam distance configuration.


Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details

Your Complementary Workflow Solution for Acquisition of Serial Sections

Up to 1000 serial sections can be collected in one day by the ATUMtome, an automated tape collecting ultramicrotome. Subsequently, the tape with sections is mounted on a silicon wafer and can be imaged with a ZEISS light microscope using ZEN imaging software and Shuttle & Find.

By taking a light microscope overview image you can plan your experiment and navigate easily on your sample within the MultiSEM using the same ZEN software user interface. For planning and setting-up the acquisition workflow only one graphical user interface is needed. Automated section detection supports you in identifying and targeting your regions of interest in a very efficient way.

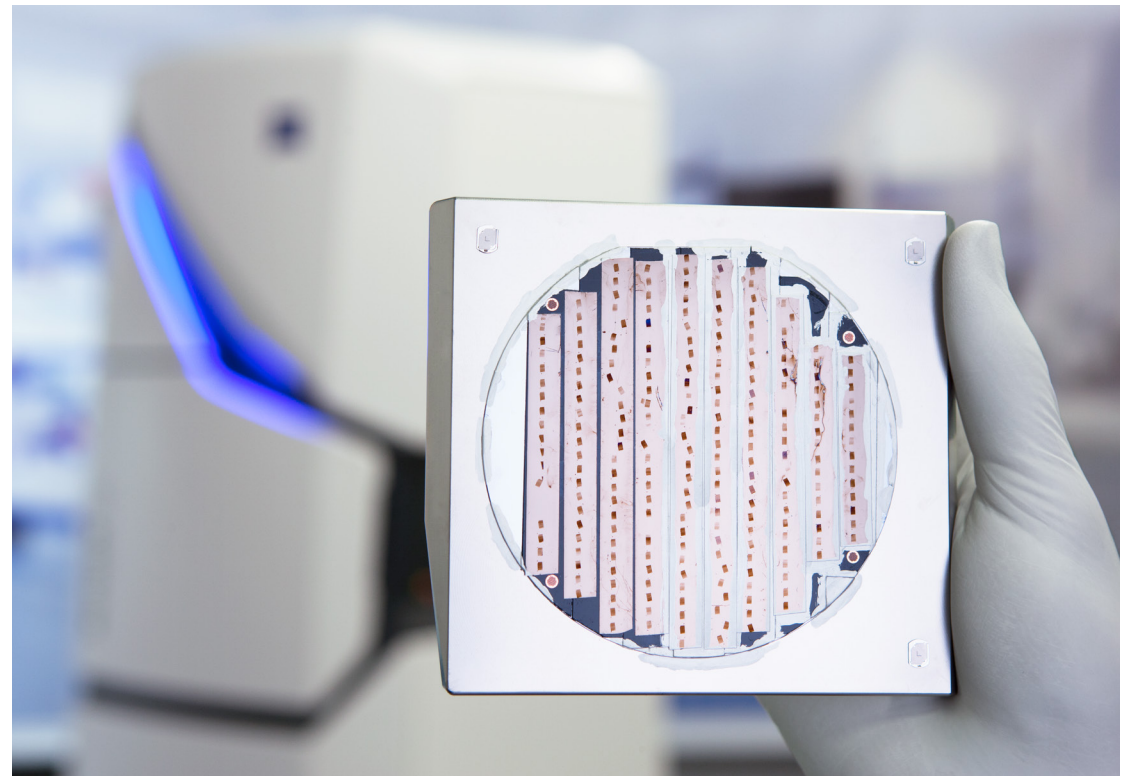

Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details

Automated Section Detection



ZEN for MultiSEM employs a section detection algorithm that saves you hours of manual labor. By template matching all sections are labeled automatically and can be defined as regions of interest for the experiment setup. Selected sub-regions can be transferred from one section to all others by just one mouse click.

Intelligent Retake Manager

ZEN for MultiSEM features a unique and intelligent data management system to ensure data completeness. Assess your data quality by using the high resolution viewer. If required, the smart retake manager supports you in setting up the follow-up experiment. Retake images are seamlessly added to the existing data set.

Image Processing and Online Analysis

The high-performance image acquisition pipeline of MultiSEM provides processing tools that run in the background during data acquisition. A flexible user interface supports image stitching for a wide range of samples.

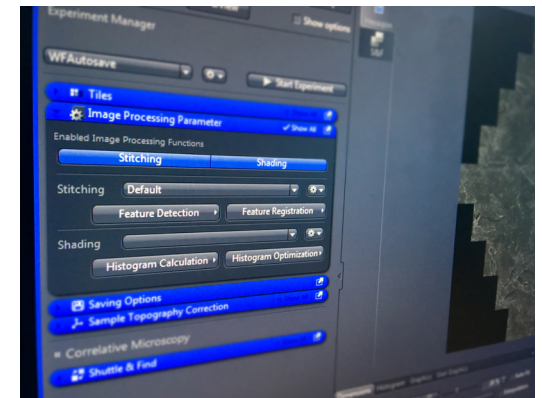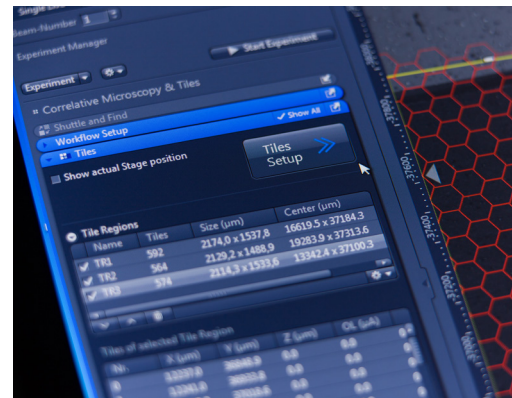


Tailored Precisely to Your Applications

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details


Typical Applications, Typical Samples	Task	ZEISS MultiSEM Offers
Stained Serial Sections of Brain Tissue	Capture large images of sections for subsequent 3D reconstruction needed for analysis in connectomics	Highest throughput electron microscopy at high resolution
Ultrathin Sections from Cultured Cells or Tissue Cultures	Screen through large sets of samples with different treatments and compare the results.	Larger regions of interest (ROI) in less time and for complete experimental trials, statistics become more reliable
Reverse engineering, Computer Chips, Patterned Silicon Wafers	Examine large areas with nanometer-sized structures	Imaging of entire chip surfaces in reasonable time frames
Analysis of Polished Rock Samples	Examine large sample surface areas to evaluate natural resources	Better quantitative assessment, larger ROIs, better statistics
Characterization of critical battery components	Analysis of large sample areas to quantify aging processes in batteries	Real, large-area quantitative assessment, delivering significant statistically relevant data
Resin-embedded bone samples	Analysis of cell distribution patterns on the etched surface of bulk bone samples	Retrieving contextual information about small features of interest at a macroscopic scale
Characterization of electronic components	Analysis of large sample areas to evaluate manufacturing processes	Real, large-area quantitative assessment, delivering significant statistically relevant data

ZEISS MultiSEM at Work


- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details

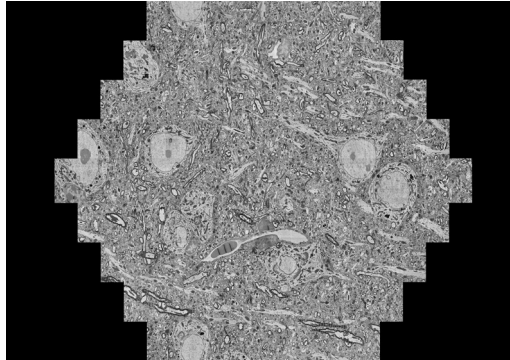


Mouse brain, 50 nm thick section. Image acquired with a MultiSEM covering a hexagonal field of view of $108\ \mu\text{m} \times 94\ \mu\text{m}$ at 4 nm pixel size. Sample: courtesy of J. Lichtman, Harvard University, Cambridge, MA, USA.

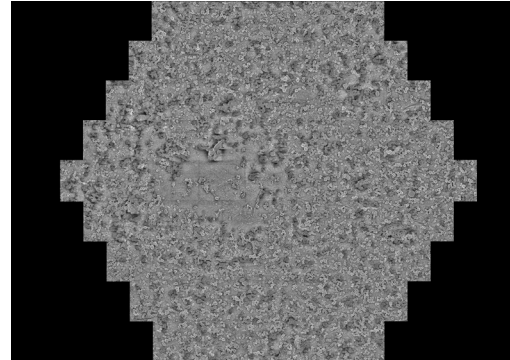


Separator foil of a cycled battery with precipitates from the anode side. Image acquired with a MultiSEM at low landing energy of 1 keV and 4 nm pixel size, covering a field of view of $108\ \mu\text{m} \times 94\ \mu\text{m}$. Sample: courtesy of U. Golla-Schindler and T. Bernthaler, Hochschule Aalen, Aalen, Germany.

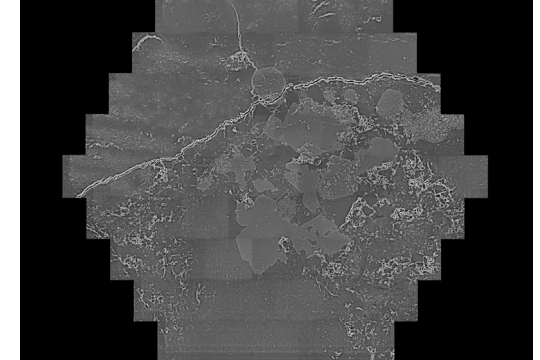


High maturity shale rock sample with broad ion beam milled surface. Image acquired with a MultiSEM at 4 nm pixel size, field of view is $108\ \mu\text{m} \times 94\ \mu\text{m}$. Sample: courtesy of L. Hathon, University of Houston, Houston, TX, USA.

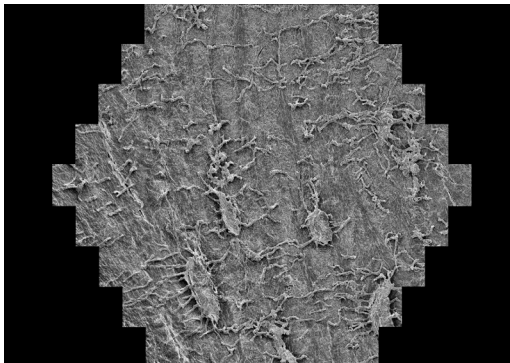

Femoral neck sample, PMMA embedded, surface-polished and selectively etched to carve out osteocytes, hidden within the bone matrix before. Image acquired with a MultiSEM at 12 nm pixel size, field of view is $135\ \mu\text{m} \times 117\ \mu\text{m}$. Sample: courtesy of M. Knothe Tate, University of New South Wales, Australia, and Ulf Knothe, Cleveland, OH, USA.

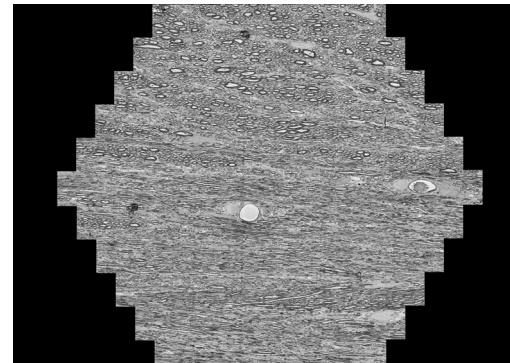

Mouse brain, 50 nm thick section. Image acquired with a MultiSEM covering a hexagonal field of view of $165\ \mu\text{m} \times 143\ \mu\text{m}$ at 4 nm pixel size. Sample: courtesy of J. Lichtman, Harvard University, Cambridge, MA, USA.

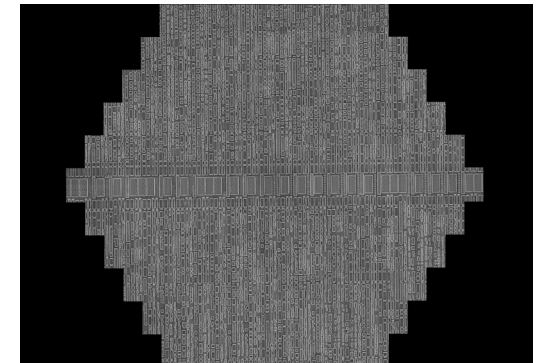

65 nm technology node graphics processor integrated circuit, stripped to its silicon substrate with HF acid etching. Image acquired with a MultiSEM at 4 nm pixel size, covering a hexagonal field of view of $165\ \mu\text{m} \times 143\ \mu\text{m}$.


Your Flexible Choice of Components

- › In Brief
- › The Advantages
- › The Applications
- › **The System**
- › Technology and Details

Accessory	Function	Technical Details
Chamber Plasma Cleaner	For cleaning the MultiSEM chamber. Reduction of contamination by hydrocarbons results in improved image quality and resolution	Generation of reactive gas-phase radicals removing unwanted contaminants. Consists of <ul style="list-style-type: none"> ■ Plasma Cleaner Evactron Zephyr ■ Adapter kit for MultiSEM chamber ■ Control integrated in ZEN software
Sample Plasma Cleaner	For cleaning and etching the sample in the airlock. Reduction of sample surface contamination results in improved image quality and resolution	Generation of reactive gas-phase radicals removing unwanted contaminants. Requires airlock. Consists of <ul style="list-style-type: none"> ■ Plasma Cleaner Evactron Zephyr ■ Multiport for MultiSEM airlock ■ Adapter Kit for Multiport ■ Control integrated in ZEN software
Standard Sample Holder	Flat surface holder for flexible sample mounting ($\leq 100 \times 100 \text{ mm}^2$)	Including L-marker fiducials for Shuttle & Find functionality
Multi-Purpose Sample Holder	For mounting standard sized EM stubs and silicon wafer chips. Additional space for flexible sample mounting (ca. $50 \times 50 \text{ mm}^2$)	Including L-marker fiducials for Shuttle & Find functionality, dedicated slots for standard sized EM stubs ($6 \times 12.7 \text{ mm}$, $3 \times 25.4 \text{ mm}$, $2 \times 32.0 \text{ mm}$), $6 \times$ silicon wafer chips ($10 \times 10 \text{ mm}$)
Multi-Purpose Sample Holder for Life Sciences	For mounting standard sized EM stubs, TEM grids and ITO cover slips. Additional space for flexible sample mounting (ca. $45 \times 35 \text{ mm}^2$)	Including L-marker fiducials for Shuttle & Find functionality, dedicated slots for standard sized EM stubs ($6 \times 12.7 \text{ mm}$, $2 \times 25.4 \text{ mm}$, $1 \times 32.0 \text{ mm}$), $6 \times$ silicon wafer chips ($10 \times 10 \text{ mm}$), $8 \times$ TEM grids and $2 \times$ cover slips
Adapter Plate for ZEISS Light Microscope	For mounting MultiSEM sample holders directly onto light microscope stage	Suitable for ZEISS Axio Imager Vario
Vibration Isolation Platform	Damping of low-frequency building vibrations (only needed if vibration specification of installation requirements is difficult to meet at customer site)	Consists of <ul style="list-style-type: none"> ■ Vibration isolation platform with active regulation by control electronics ■ Modified MultiSEM housing which enables microscope installation on vibration isolation platform
Workflow Add-ons	Function	Technical Details
Light Microscope	Large area imaging ($> 10 \text{ cm}^2$) for fast sample overview and region of interest selection. Sample positions can be accurately identified and relocated within MultiSEM	Recommended microscope is ZEISS Axio Imager.A2 Vario. ZEN software (blue edition) with ZEN connect licence required.
ATUMtome	Automated sectioning and section collection of resin-embedded biological tissue	Ultramicrotome based section collection robot from RMC Boeckeler. Up to 1000 sections per day with a typical sample thickness of $30 - 50 \text{ nm}$

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**

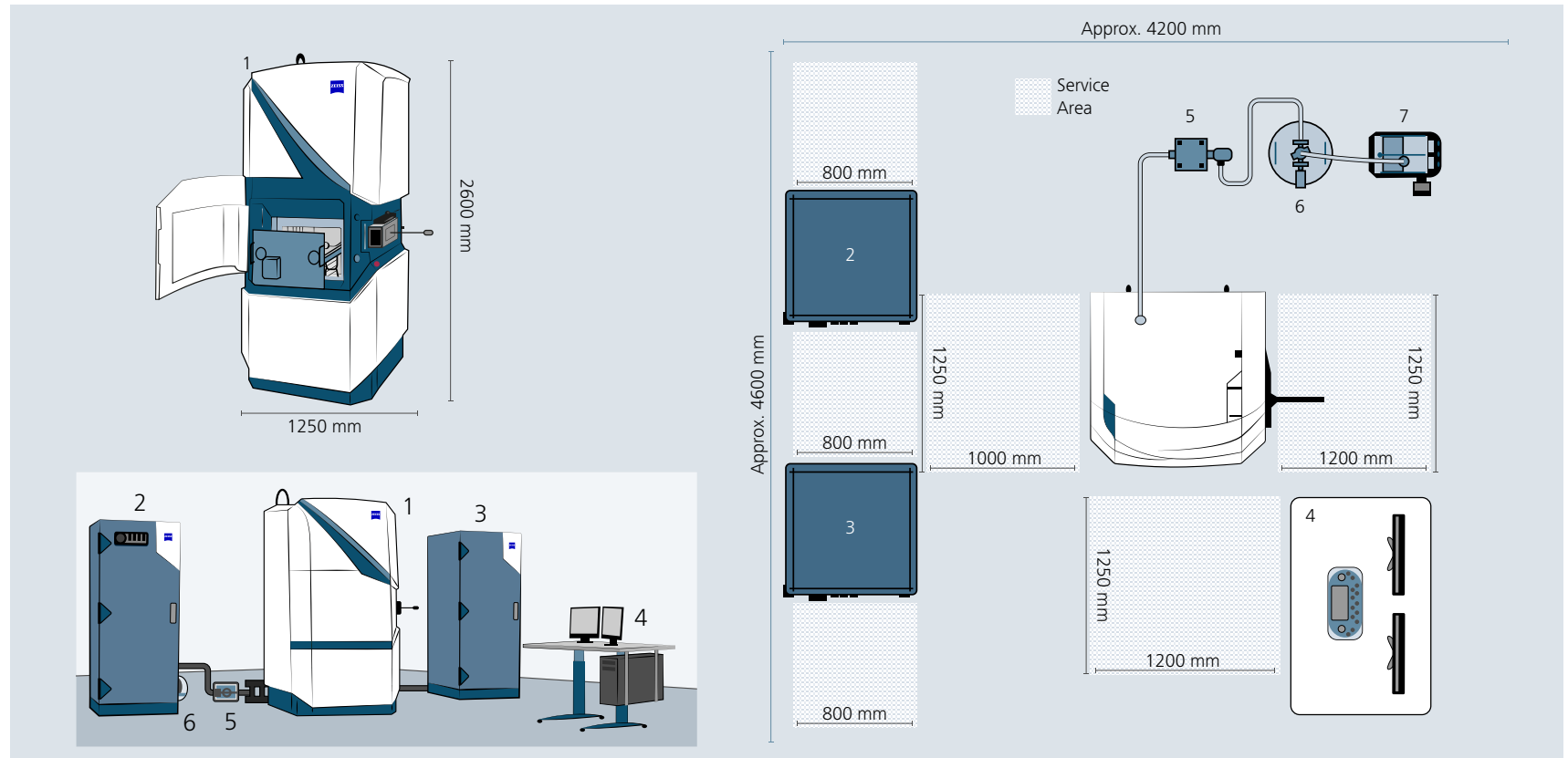

No.	Description	Size (mm) approx.	Weight (kg) approx.	Load Distribution (kg) approx.	Footprints (mm) approx.
1	Microscope	1250 × 1200 × 2600	1600	4 × 400	4 × ø 156
2	PSU Rack	800 × 800 × 2100	300	4 × 75	4 × ø 39
3	IA Rack	800 × 800 × 2100	332	4 × 83	4 × ø 39
4	Table and PC	1153 × 980 × 1350	112	4 × 28	4 × ø 50
5	Static Damping Block	180 × 180 × 160	13	1 × 13	180 × 180
6	Quiet mode	350 × 350 × 550	5	1 × 5	350 × 350
7	Pre-Vacuum Pump	430 × 250 × 290	25	1 × 25	200 × 180

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**

Electron Optics		MultiSEM 506
Resolution	Average resolution of all beams @ 1.0 kV, 1.5 kV, 3.0 kV	≤ 3.5 nm
Resolution Uniformity	@ 1.0 kV, 1.5 kV, 3.0 kV and 12 μ m pitch size	$\leq \pm 1$ nm
Landing Energy	Range	1.0 – 3.0 kV
Beam Arrangement	Beam pattern	Hexagonal
	Number of beams	91
	Pitch size (width of single beam image)	12 μ m or 15 μ m (optional)
	Pitch uniformity	$\leq \pm 1$ %
Field of View (FoV)	Long axis of hexagon (12 μ m pitch / 15 μ m pitch)	132 μ m / 165 μ m
Beam Current	Single beam	≥ 570 pA
	Total current	≥ 52 nA
	Uniformity	$\leq \pm 10$ %
Electron Source	Filament	Schottky emitter
	Filament current stability	≤ 1 % per hour
Beam Blanker		Electrostatic beam blanker
Working Distance		1.4 mm
Detection		Secondary electron projection optics with high-efficiency multi detection unit
Scanning		
Scan Rate		Max. 20 MHz per beam, different discrete scan speeds are available
Scan Mode		Step and scan
Pixel Size	Range at 12 μ m pitch	2 nm – 20 nm
	Range at 15 μ m pitch	2.5 nm – 20 nm
Scan Arrangement		Image tile consists of 91 Sub-Images arranged in a hexagonal pattern
		Adjustable overlap of adjacent scan fields

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**

Stage and Specimen		MultiSEM 506
Stage	Type:	Stepper Stage
	Usable travel range x/y/z:	100 / 100 / 30 mm
	Repeatability XY	$\leq \pm 3 \mu\text{m}$
	Settling time	$\leq 1.5 \text{ s}$
Specimen Requirements	Maximum size in XY	$100 \times 100 \text{ mm}^2$
	Maximum height	$\leq 30 \text{ mm}$
	Maximum flatness	$\leq 500 \text{ nm} / 100 \mu\text{m}$ (Peak-to-Peak)
	Maximum weight	$\leq 0.2 \text{ kg}$
Specimen Exchange Time	With airlock	$\leq 5 \text{ min}$
Software		
User Interface		ZEN for MultiSEM
Application Programming Interface (API)		Provided for custom workflow development
Shuttle & Find Functionality		Reliable transfer of sample coordinates from different imaging modalities (e.g. light microscope or single-beam SEM)
Performance Monitoring		Check of all relevant system parameters such as vacuum pressures or optics alignment quality
Parallel Software Architecture		Distributed image acquisition
Data Base Support		Provided for workflow and data management
Automated Alignment Functions		Autofocus, auto-stigmation, detector equalization, etc.
Image Acquisition Workflow		
Graphical Experiment Setup		Image based region of interest selection
Automated Section Detection (Option)		Fast workflow setup for serial sections imaging
Interaction Requirement		Max. 1 hr/24 hrs dedicated user interaction for beam alignment & calibration

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**

Computer Hardware		MultiSEM 506
Main Controller		≥ 4 core CPU (64 bit), ≥ 32 GB DDR, ≥ 1 TB HD, min. 2 ports with 1 Gbit Ethernet
Main Controller Operating System		Windows® 10 (64 bit)
Display		2 Monitors, 1920×1080 Pixel, 24"
Image Acquisition		12 PCs, in 19" Rack
Image Acquisition PC		≥ 4 core CPU (64 bit), ≥ 32 GB DDR3, 1 Gbit Ethernet, 10 Gbit Ethernet to customer network
Data Transfer Rate to Customer Network		≥ 10 Gbit Ethernet
Vacuum System		
Chamber Vacuum Pumps		Turbo molecular pump (≥ 260 l/s); oil-free Scroll Pump
Chamber Operating Pressure		≤ 1 × 10 ⁻⁵ mbar
Monitoring		Automatic monitoring of all relevant pressures



Carl Zeiss Microscopy GmbH
07745 Jena, Germany
microscopy@zeiss.com
www.zeiss.com/multisem

Follow us on social media:

