
EyeMag Pro S

Instructions for Use
Εγχειρίδιο λειτουργίας
Használati utasítás
Instrukcja obsługi
Instrucţiuni de utilizare
Návod na používanie
Navodila za uporabo
Kullanım Kılavuzu

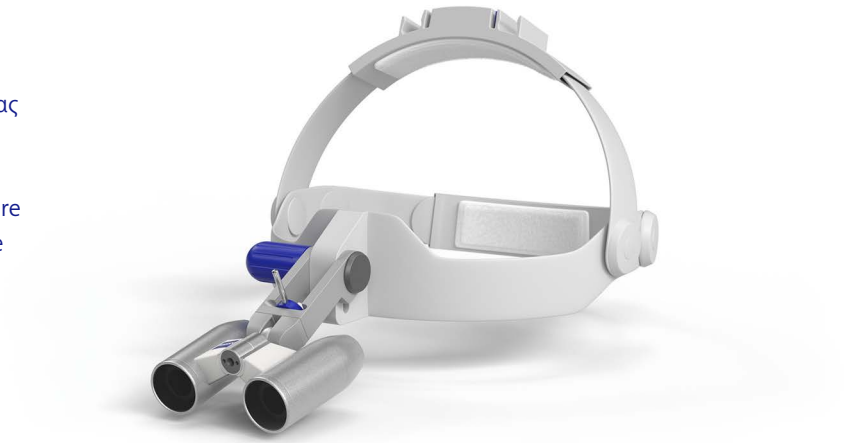Congratulations on purchasing your medical head-worn loupe from ZEISS. To ensure continuous and trouble-free use, please observe the information provided in these Instructions for Use.

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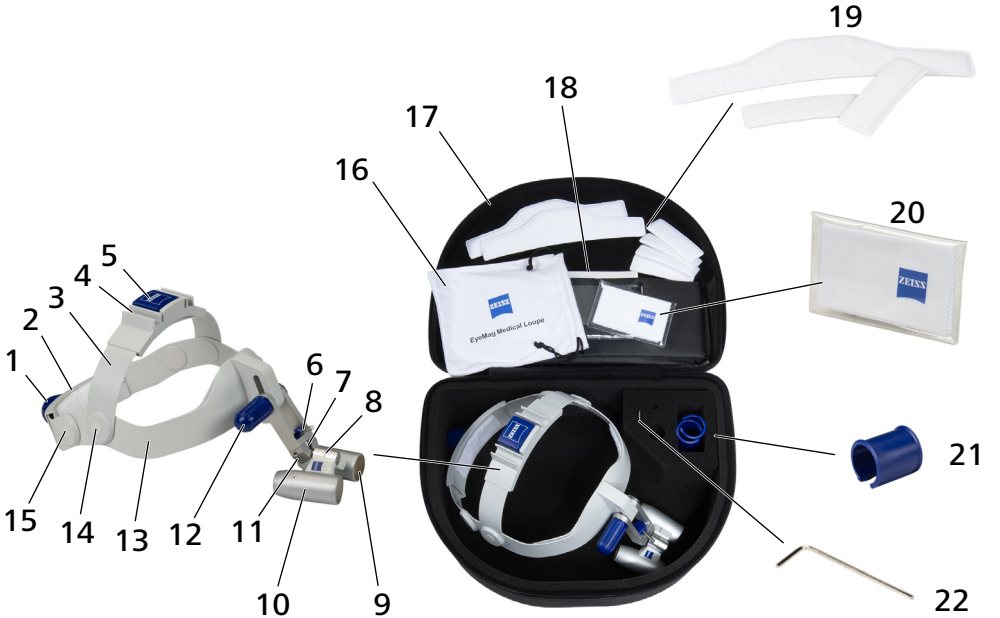
EN

Table of Contents

Scope of supply 5
Safety measures 5
Hazard symbols 5
Product design fundamentals 6
Regulatory information 6
Target group 6
Intended use 6
Typical misuse 6
Notes for the user 6
Prevention of injury 7
Prevention of damage to material 7
Other safety measures 7
Notifying the manufacturer and authorities 7
Figures 7
Symbols and labels 8
Description 9
Product 9
Components 9
Operation 10
First steps 10
Configuring the product 11
Configuring the optical system 13
Flip-up function 15
Replacing the optical system 16
Replacing the textile pads 17
Using components 17
Protection caps 17
Contact guard 18
Using the product with LED illumination 19
Care and maintenance 19
Storage 19
Cleaning 19
Disinfection 20
Disposal 20
Product data 20
Instructions for Use 21
Accessories 21
ZEISS LED illumination 21
System carrier 21
Optical systems 21
Technical specifications 22

EN

Scope of supply

The following parts belong to your EyeMag Pro S (hereinafter referred to as "product"):

1 Head circumference adjustment knob
2 Rear headband (with attached textile pad)
3 Top headband (with attached textile pad)
4 Notch for clamping sleeve (2x)
5 Adjustment knob for top headband
6 Locking lever
7 Ball adapter
8 Folding bridge
9 Lens protection caps (mounted)
10 Optical system
11 Optical mount
12 Height adjustment locking knob
13 Front headband (with attached textile pad)
14 Hinge for top headband
15 Hinge for rear headband
16 Fabric bag
17 Softcase, large
18 Instructions for Use G-30-2331 & G-30-1560
19 Additional replaceable textile pads
20 ZEISS Microfiber Lens Cleaning Cloth
21 Contact guard (2x)
22 SW1.5 hex key

Technical details and delivery package subject to change.

Safety measures

Correct use of the device is absolutely vital for safe operation. Before using the product, please note the information contained in these Instructions for Use, the Instructions for Use of any additional product equipment such as an illumination system, and the legal accident prevention regulations. These Instructions for Use are part of the product and must be accessible at all times. Keep the Instructions for Use over the entire service life of the product.

Hazard symbols

The following hazard symbols and signal words have been incorporated in the safety information of these Instructions for Use.

Symbol	Signal word	In the event of noncompliance
⚠	**WARNING!**	**Death** or **severe injuries** may occur
⚠	**CAUTION!**	**Minor** or **moderate** injuries may occur
	NOTE!	**Damage to property** may occur

EN

Product design fundamentals

The product described in these Instructions for Use has been designed and tested in accordance with ZEISS safety standards as well as national and international regulations.

The product as well as the components are latex-free.

The product consists of an optical system (Keplerian binocular teleloupe system) and a system carrier.

Carl Zeiss Meditec AG uses an EN ISO 13485 compliant QM system certified by DQS.

Regulatory information

In the medical field, the loupe is a Class 1 product in accordance with EU Medical Device Regulation (MDR) 2017/745.

Target group

These Instructions for Use are intended for medical personnel responsible for the use, cleaning and disposal of the product.

Intended use

The loupes of the EyeMag family are medical loupes used for binocular stereoscopic near vision, e.g., in dentistry or surgical applications.

Typical misuse

Do not use the product

– for looking into bright light sources, e.g. lasers
– if modifications have been made to the product

Notes for the user

The following information symbols are used in these Instructions for Use:

Symbol	Application
–	Bulleted list
1.	Call to action where the order of steps must be followed.
■	Prompt for an action

The product causes hazard when not correctly used or when not used in a proper condition. Please note the information provided below.

Prevention of injury

WARNING!

Intense light radiation may lead to irreversible damage to the eye!

- When wearing the product, never look directly at the sun or at powerful technical light sources such as lasers.
- Use the product only for the purpose described.

Caution!

Risk of crushing fingers due to incorrect operation!

- Please note the steps of operation described in these Instructions for Use.

Falling parts, body fluids or contamination can infect patients or users.

- Before use, check the product for damage, loose parts and for proper function.
 - Do not continue using the product if it is damaged.
 - The parts that you can replace yourself can be found under "Accessories", on page 21.

Prevention of damage to material

Mechanical stress will damage the product!

- Do not subject the product to severe mechanical stress, such as by dropping it or knocking it against other objects.

The ingress of liquid will damage the product!

- Do not clean the product in an ultrasonic bath, under running water, in a water bath or in a similar manner, as it is not waterproof.

Other safety measures

- If it is absolutely necessary to use the product for a surgical procedure, perform the procedure only if a substitute product of the same type is available.

Notifying the manufacturer and authorities

- If a serious incident occurs in connection with this medical device affecting the operator or another person, the operator (or responsible person) must report this serious incident to the manufacturer or distributor of the medical product. In the European Union, the operator must report this serious incident to the responsible authority in the country in question.

Figures

Figures may differ from the original.

EN

Symbols and labels

Symbol	Description	
	**UDI label** – Manufacturer (company name) – Symbol for "Manufacturer" – Manufacturer's address – Product name – Machine readable label (QR code) – Reference number – Serial number – Date of manufacture – MD designates the product as a medical device – UDI mark – UDI Device Identifier – UDI Production Identifier – CE mark	 Carl Zeiss Meditec AG Goeschwitzer Strasse 51-52 07745 Jena, Germany EyeMag Pro S REF SN MD UDI UDI-DI UDI-PI CE
X,Xx / XXX	Magnification factor / working distance in mm	
ZEISS	ZEISS logo (manufacturer)	

Fig. 1

Fig. 2

Fig. 3

Description

Product

The product (Fig. 1) consists of an S system carrier (support) and a binocular optical system. The system carrier is comprised of flexible hinges and straps. The optical system features two eyepieces. The optical system provides large homogeneous fields of view, giving you an optimal overview of the treatment area. Information regarding the magnification and working distance of the optical system is provided on the bridge between the sleeves.

Components

The following components can be attached to the product:

Protection caps (Fig. 2)

The protection caps protect the optical system against contamination and mechanical damage during use.

The product is shipped with the protection caps already attached to the optical system.

Contact guard (Fig. 3)

The sterilizable contact guard permits sterile operation of the flip-up function.

EN

Operation

The product allows the user to work at all times with his or her head and body in a relaxed and comfortable upright position.

The various configuration and attachment options available are described on the pages to follow.

First steps

WARNING!

Risk of infection due to loose or falling parts of the product!

- Before use, check the product and its components for damage, loose parts and functional safety.

- Adjust the product's settings as described in this manual.

Fig. 4

Configuring the product

To obtain optimum wearing comfort, you must adjust the product to fit the size and shape of your head.

The product generally only needs to be configured once. For safety reasons, this should be done in front of a mirror.

Proceed as follows:

1. Flip the optical system up into the last clickstop position.
2. Tighten the locking knob so that the optical system remains in the configured position.
3. Turn the adjustment knob to increase the head circumference until the product can be comfortably placed on your head.
4. Lengthen the top headband until the product can be comfortably placed on your head.
5. Put on the product.
6. Lower or raise the rear headband as necessary.
7. Rotate the top headband into the desired clickstop position as necessary.
8. Turn the adjustment knob to decrease the head circumference until the product sits firmly on your head.
9. Align the product so that the optical mount sits symmetrically above the bridge of the nose, and the lower edge of the front headband is approx. 35 to 45 mm above the center of the eyes.

EN

10. Push the two pieces of the top headband together to reduce the head height.
11. To increase the head height, press the head height adjustment button and pull the two pieces of the upper headband apart.
12. To remove the product, turn the head circumference adjustment knob and lift the product upwards.

Fig. 5


Fig. 6


Configuring the optical system


If you have adjusted the product to fit the size and shape of your head as described on the previous pages, the optical system should be in the upper position.

1. Flip the optical system into the lower position.
 When doing this, ensure that there is sufficient space between the system and your eyes, and avoid direct contact between the optical mount and the bridge of your nose.
2. Hold the optical system firmly with one hand.
3. Turn the height adjustment locking knob to the left to loosen the height adjustment lock.
4. Move the clamping lever all the way to the left to release it.

EN

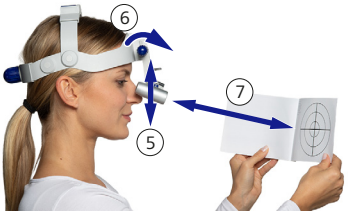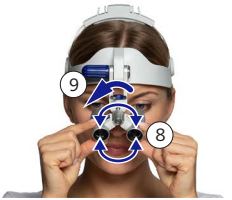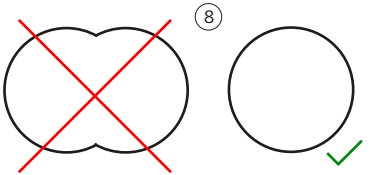5. Adjust the optical system to your desired height.
6. Turn the height adjustment locking knob to lock the height position.
7. Adjust the distance to the test card (Fig. 16) until the test card is in the focal plane of the optical system.
8. Adjust the foldable bridge to set the interpupillary distance. If you see a circular field of vision in the focal plane, you have found the optimal setting.
9. To lock your configuration settings, move the clamping lever all the way to the right.

Fig. 7

Fig. 8

Fig. 9

Flip-up function

If you do not need to use the product while wearing it, the optical system can be flipped up and out of the visual field (Fig. 8).

NOTE!

Quick and jerky movements may damage the carrier, optical system or eyeglass lens!

- Always be careful when flipping the optical system up or down (Fig. 8, 9).

If you are using a mounted sterilizable contact guard, please observe the following when using the flip-up function:

WARNING!

Risk of infection caused by non-sterile surfaces!

- Do not touch any unsterile surfaces. If you can no longer ensure that the contact guard is sterile, have another person operate the flip-up function.

EN

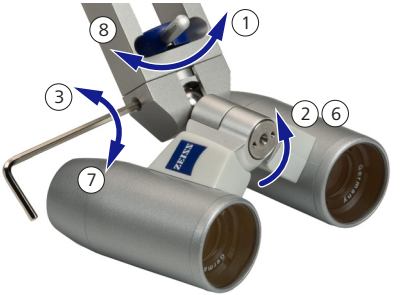
Fig. 10

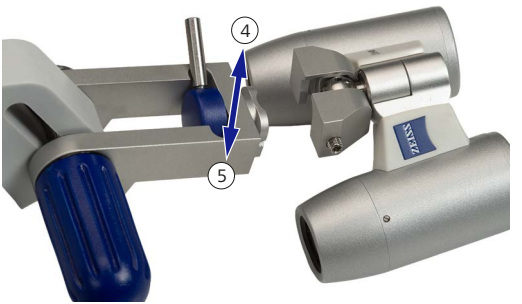
Fig. 11

Replacing the optical system

(Fig. 10 and 11)

WARNING!

Risk of infection due to loose or falling parts of the product!

- Before use, check the product and its components for damage, loose parts and functional safety.
- If you are unable to restore functionality, do not continue to use the product.

NOTE!

Risk of damage to the optical system or system carrier caused by falling!

- Firmly hold the product with both hands when removing the optical system.

1. Move the clamping lever all the way to the left to release it.
2. If you need to, flip the optical system out of the way.
3. Use the SW1.5 hex key to loosen the set screw located on the side of the optical mount.
4. Push the optical mount and optical system upward along the guideway to remove them.
5. Insert an optical system of your choice along the guideway from above to attach it to the optical mount.
6. If you need to, flip the optical system out of the way.

7. Use the SW1.5 hex key to tighten the set screw located on the side of the optical mount.
8. To lock the optical system in place, move the clamping lever all the way to the right.
9. Configure the optical system as described under "Configuring the optical system" on page 13.

Replacing the textile pads

The support is outfitted with textile pads in various positions to improve wearing comfort and hygiene (Fig. 12). These replaceable textile pads are attached by means of velcro fasteners. They are washable (delicate cycle at 40 °C).

Fig. 12

Using components

Protection caps

Removing the protection caps (Fig. 13)

- Carefully unscrew the protection caps from the optical system.

Attaching the protection caps (Fig. 13)

- Place the protection caps on the optical system and carefully screw them on.

Fig. 13

EN

Fig. 14

Fig. 15

Contact guard

Attaching the contact guard

WARNING!

Risk of infection caused by unsterile contact guard!

- The sterilizable contact guard permits sterile operation of the flip-up function.
- Do not touch any unsterile surfaces. If you can no longer ensure that the contact guard is sterile, have another person operate the flip-up function.
- The contact guard must be sterilized before each use. This also applies to the first use after delivery. See G-30-1560.

- Slide the already sterilized contact guard over the front of the optical system (Fig. 14).

Removing the contact guard

- Slide the contact guard towards the front of the optical system to remove it (Fig. 15).

Using the product with LED illumination

The product can be used in combination with ZEISS LED illumination (see page 21). A corresponding adapter is supplied with the purchase of the LED illumination, or can be ordered separately.

Care and maintenance

The product (including its components) is a sensitive optical product and must be handled with care. Therefore, please note the following recommendations.

Storage

Store the product in its original case (large softcase) when it is not in use. Make sure that no parts come under pressure.

Cleaning

WARNING!

Risk of infection due to surface parts coming loose and falling off!

- Do not use solvents or sharp and grinding cleaning agents as this may cause particles to be released from the surface.

CAUTION!

Incompatibility of cleaning agents!

- To avoid allergic reactions, follow the instructions provided by the manufacturer of the cleaning agent.

NOTE!

The product is not waterproof! Penetrating liquid may damage the product.

- Do not clean the product in an ultrasonic bath.
- Do not clean the product under running water.
- Do not clean the product in a water bath.
- Do not spray cleaning fluids directly onto the product.

- Clean the product with standard cleaning agents after every use.
- Clean the contaminated glass surfaces using the cleaning towel included in the delivery package.

EN

Disinfection

NOTE!

Surface damage to the product!

- Use an aldehyde and/or alcohol-based disinfectant. The addition of quaternary compounds is acceptable. To avoid damaging the surfaces, do not use any disinfecting components other than those specified below.

The maximum concentrations are:

– For alcohol (tested with 2 propanol): 80%
– For aldehyde (tested with glutaraldehyde): 2%
– For quaternary compounds (tested with DDAC): 2.5%
– For phenols (tested with phenylphenol): 0.7%

Disposal

NOTE!

Incorrect disposal contaminates the environment!

- Observe the local laws/regulations governing disposal.

- Keep the packaging material for use in the event of repair.
- If you want to dispose of the packaging materials: Dispose of the packaging materials by sending them to an approved collection system for recycling.

Product data

Use the product only with the components and accessories included in the delivery package and approved by ZEISS. The components and accessories listed below can be ordered from ZEISS. You can find the ZEISS representative for your country online at:
www.zeiss.com/med/contacts

Instructions for Use

Description	Order number	Pcs.
EyeMag Pro S G-30-2331	000000-2801-502	1
Reprocessing Instructions G-30-1560	000000-1312-845	1

Accessories

Description	Order number	Pcs.
Contact guard	304156-9900-000	2
ZEISS Lens Microfiber Cloth	000000-0576-756	1
Softcase, large	000000-1460-441	1
Fabric bag	304100-9001-000	1
Textile pads, rear and upper	304160-8011-000	1
Textile pad, front	304160-8012-000	1
SW1.5 hex key	000000-0015-244	1

ZEISS LED illumination

Description	Order number	Pcs.
EyeMag Light II	304121-9010-000	1

System carrier

Description	Order number	Pcs.
System carrier, S (support)	304160-9100-000	1

Optical systems

Description	Order number	Pcs.
K 3.2 x/ 500 mm	304156-9503-000	1
K 3.3 x/ 450 mm	304156-9453-000	1
K 3.5 x/ 400 mm	304156-9403-000	1
K 3.6 x/ 350 mm	304156-9353-000	1
K 4.0 x/ 300 mm	304156-9304-000	1
K 4.0 x/ 450 mm	304156-9454-000	1
K 4.0 x/ 500 mm	304156-9504-000	1
K 4.0 x/ 550 mm	304156-9554-000	1
K 4.3 x/ 400 mm	304156-9404-000	1
K 4.5 x/ 350 mm	304156-9354-000	1
K 5.0 x/ 300 mm	304156-9305-000	1
K 5.0 x/ 500 mm	304156-9505-000	1
K 5.0 x/ 550 mm	304156-9555-000	1
K 6.0 x/ 500 mm	304156-9506-000	1
K 6.0 x/ 550 mm	304156-9556-000	1
Protection caps	304156-9901-000	2

EN

Technical specifications

Working distance in mm (measured from the eye)	500	450	400	350	300	450	500	550	400	350	300	500	550	500	550
Optical system magnification	3.2x	3.3x	3.5x	3.6x	4.0x	4.0x	4.0x	4.0x	4.3x	4.5x	5.0x	5.0x	5.0x	6.0x	6.0
Approx. max weight of optical system	93 g														
Approx. max weight with STMS system carrier	330 g														
Average service life	6 years														

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