



## VISUSCOUT 100 from ZEISS

[www.zeiss.com/visuscout100](http://www.zeiss.com/visuscout100)



Seeing beyond

# Fundus imaging – wherever you need it

## Truly portable imaging from ZEISS

Reliably detecting and monitoring retinal disorders is key to ensuring high-quality care and to maintaining the vision of your patients. The VISUSCOUT® 100 from ZEISS lets you do precisely that. As a mobile fundus camera, it is the perfect imaging companion.



## Retinal diagnosis that reaches every patient

With its lightweight and compact design, the ZEISS VISUSCOUT 100 allows simple retinal imaging and great flexibility when and wherever the need arises.

- For checkups and documentation of retinal conditions such as diabetic retinopathy, glaucoma and AMD
- As a secondary device for immobile patients
- Ideal for outreach programs, mobile eye care and satellite offices
- When working with children
- For general consultations and home visits
- As a simple fundus camera for small practices with space constraints
- Well-suited for health care professionals interested in an easy, cost-effective fundus imaging solution

Packed into a small, rugged carrying case, the ZEISS VISUSCOUT 100 can be conveniently transported and easily fits into any practice setup. Thanks to the camera's non-mydratric operation and precise autofocus function, dilation of the eyes is not required. Its battery power provides added flexibility. The WiFi functionality enables instant transfer of images to a PC.

**The ZEISS VISUSCOUT 100 comes with an ergonomic charging station. A slit lamp adapter is available as an option.**

## Highlights at a glance

### Image quality

Take color and red-free images instantly with the 40° field of view. The device fulfills all relevant ISO 10940 fundus camera standard requirements.

### Easy capture

With its on-screen targeting aid and integrated autofocus function, using the device is easy and convenient. Nine internal fixation LEDs help align the patient correctly and also facilitate the capture of peripheral images. Thanks to its battery-powered, lightweight design, system handling is not restricted by cables.

### Built-in simplicity


Operating the ZEISS VISUSCOUT 100 is simple and straightforward following a brief training session. Examinations are quickly performed and reviewed without needing to move the patient. Even medical staff can operate the device, enabling cost-effective image capture prior to the examination.

### Included PC-software

The ZEISS VISUSCOUT 100 uses standard file formats that allow you to view images and videos with widely available, off-the-shelf software solutions. When combined with the supplied PC-software, the device can be connected to a DICOM-based patient management system such as FORUM® from ZEISS. FORUM allows you to access patient work lists. It also enables you to transfer the data to a central archive.

## Key benefits for you

- Wireless flexibility for fundus imaging wherever you need it
- Easy capture of images and videos
- Integrates with DICOM-based data management systems using supplied PC-Software
- Simple to operate, even for medical staff and assistants

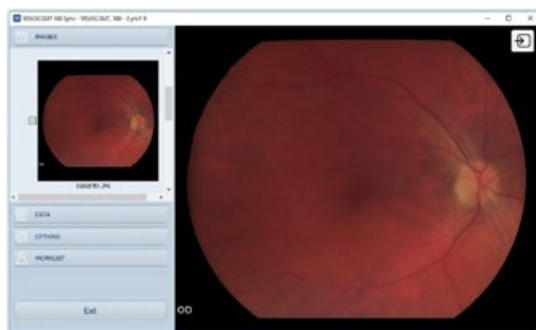

# Specifications and technical data

## VISUSCOUT 100 from ZEISS


<b>Dimensions (W x D x H)</b>	115 x 216 x 199 mm
<b>Weight</b>	800 g
<b>Sensor resolution / type</b>	5 MP / CMOS
<b>Protection class</b>	1
<b>Charging station</b>	100 - 240 V ~0.4 A 47 - 63 Hz, Output: 9 V, 1.1 A, 10 W
<b>Battery</b>	Rechargeable Li-Ion 3.6 V, 2.350 mAh
<b>Field of view</b>	40°
<b>Capture mode</b>	Color, red-free, IR
<b>Focus range</b>	-20 D to + 20 D
<b>Fixation</b>	9 x LED, internal
<b>Minimum pupil size</b>	3.5 mm
<b>Display</b>	4", TFT-LCD, 800 x 480 px, 16.7 M colors, anti-glare coating
<b>Data connectivity</b>	USB WiFi (802.11 b/g/n) on 2.4 GHz band
<b>Data formats</b>	JPEG, MPEG4/1 (DICOM)
<b>Memory</b>	WIFI SDHC card



*The ZEISS VISUSCOUT 100:  
flexible and mobile fundus imaging*

CE 0598

**Carl Zeiss Meditec AG**  
Goeschwitzer Strasse 51–52  
07745 Jena  
Germany  
[www.zeiss.com/visuscout100](http://www.zeiss.com/visuscout100)  
[www.zeiss.com/med/contacts](http://www.zeiss.com/med/contacts)

 **OPTOMED Oyj**  
Yrtytipellontie 1  
FI-90230 Oulu  
Finland

**EN\_33\_022\_0010V** International edition: Only for sale in selected countries.  
The contents of the brochure may differ from the current status of approval of the product or service offering in your country. Please contact our regional representatives for more information.  
Subject to changes in design and scope of delivery and due to ongoing technical development. VISUSCOUT and FORUM are either trademarks or registered trademarks of Carl Zeiss Meditec AG or other companies of the ZEISS Group in Germany and/or other countries.  
© Carl Zeiss Meditec AG, 2019. All rights reserved.